


1. Chairside



2. Cleaning

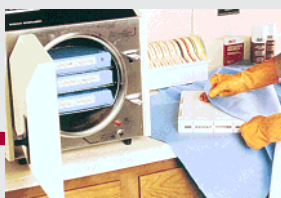


3. Sterilization



4. Storage

STAINLESS STEEL CASSETTES

IMS™ *Signature Series*™-

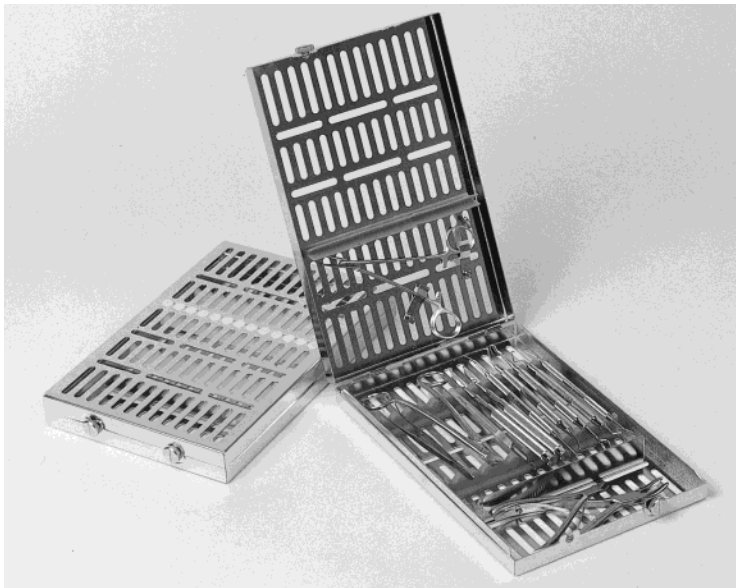

DIN Cassette for 13 Instruments + accessory area


measurements: 290 x 184 x 34 mm, mit space for e.g. 2 Langenbeck Retractors. (The mentioned Instruments belong to the Kramer-Nevis-Set).

IMDIN13

- width choice:
3 sizes
6 fittings
5 colours
- Stainless Steel for easy maintenance
- Inner space for own organisation
- Coloured Silicone rails facilitate the organisation

DIN-Cassette for 13 Instruments
5 different colour codes available

Grey	552518	IMDIN131
White	210617	IMDIN132
Orange	554758	IMDIN133
Yellow	792535	IMDIN135
Blue	130427	IMDIN138

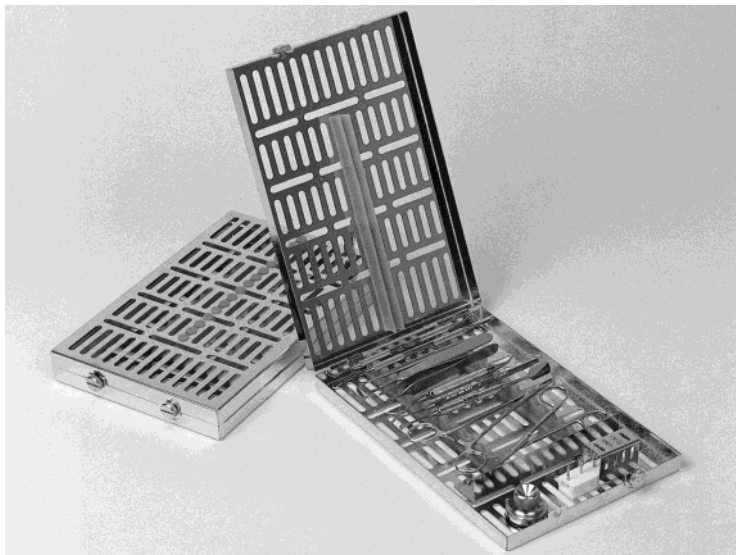

DIN Cassette for 16 Instruments + accessory area

measurements: 290 x 184 x 34 mm for Amalgam, Composite, Endodontic, Prothetik, Prophylaxe The illustration shows instruments from the Labanca-Set.

IMDIN16

DIN Cassette for 20 Instruments

measurements: 290 x 184 x 34 mm for Surgical (no illustration)

IMDIN20

DIN-Cassette for 16 Instruments
5 different colour codes available

Grey	691453	IMDIN161
White	907154	IMDIN162
Orange	280305	IMDIN163
Yellow	296899	IMDIN165
Blew	328315	IMDIN168

DIN-Cassette for 20 Instruments
5 different colour codes available

Grey	945544	IMDIN201
White	969818	IMDIN202
Orange	154571	IMDIN203
Yellow	936316	IMDIN205
Blew	496472	IMDIN208

In ordering please mention order number and colour code (last number e.g. **IMDIN208**) (see table)

- IMS Cassettes
- Cassettes IMS
- Cassette IMS
- Estuches IMS

WÜBTEN SIE ...

... that the cassettes are available in three sizes?
In DIN-size, 1/2-DIN-size and in the special measurements for STATIM-autoclaves.

SIGNATURE SERIES CASSETTES



1/2 DIN
Cassette
for 8 Instruments
+ accessory
area

IM12DIN8

measurements: 152 x 184 x
34 mm
for Diagnostic and Prophylaxe.

1/2 DIN
Cassette
for 10
Instruments

IM12DIN10


measurements: 152 x 184 x
34 mm
(no illustration)

1/2 DIN-Cassette for 8 Instruments
5 different colour codes available

Grey	746754	IM12DIN81
White	859214	IM12DIN82
Orange	136204	IM12DIN83
Yellow	845795	IM12DIN85
Blue	898252	IM12DIN88

1/2 DIN-Cassette for 10 Instruments
5 different colour codes available

Grey	358673	IM12DIN101
White	215724	IM12DIN102
Orange	858684	IM12DIN103
Yellow	266604	IM12DIN105
Blew	526142	IM12DIN108

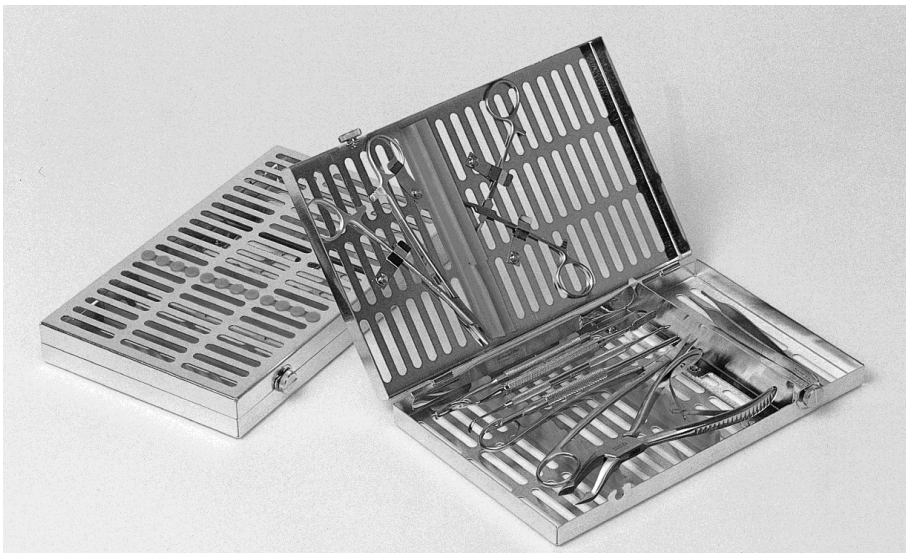

Statim Cassette
for 12
Instruments

IM812

measurements: 268 x 176 x 34 mm,
for STATIM-autoclaves.
The illustration shows instruments from
the Hourigan-Set.

Special size for STATIM-autoclaves,
for 12 Instruments,
5 different colour codes available

Grey	284467	IM8121E
White	966606	IM8122E
Orange	480975	IM8123E
Yellow	732478	IM8125E
Blew	901555	IM8128E

In ordering please mention order number and colour
code(last number e.g. IMDIN208) (see table)


- IMS Cassettes
- Cassettes IMS
- Cassette IMS
- Estuches IMS

PFLEGETIP

Our Cassettes are made from best quality
Stainless Steel and withstand
heat to 176°C.

ACCESSORIES

IMS™ NEEDLE CAP HOLDER

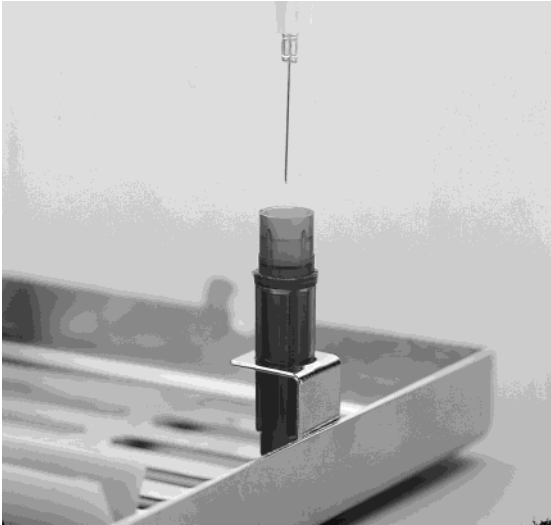


IMS™
2 Needle
Cap Holder
719145
IM1001

IMS™
6 Needle
Cap Holder
170224
IM1003

The needle cap holder allows for efficient one-handed needle recapping. Fits all conventional cover sizes. Maintains needle cap in upright position for optimum chairside efficiency. The "snap-in" assembly allows for easy placement throughout the entire cassette. Can be easily inserted into any Signature Series Cassette. (Set of 6 Needle Cap Holder + 6 Instrument Clips = **IM1004** - 220558)

IMS™ INSTRUMENT CLIP

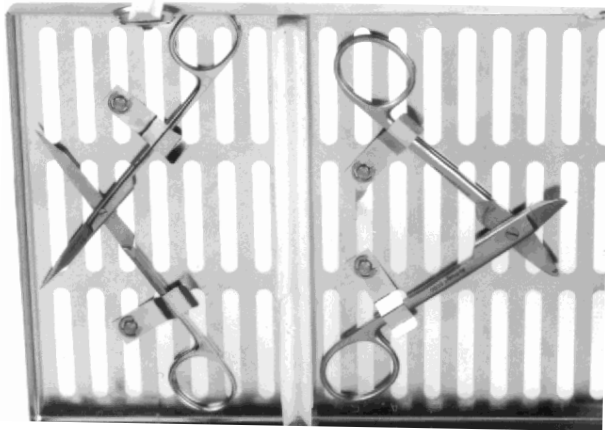

IMS™ 2 Clips
for hinged
instruments
776491
IM1000

IMS™ 6 Clips
for hinged
instruments
505358
IM1002

Hinged instrument clips secure scissors, hemostats and qir/water syringe tips for maximum efficiency and set-up organization. Hinged instruments are securely held in open position allowing for proper cleaning, sterilization and optimum instrument organization. Fits on any Signature Series Cassette.

IMS Accessories

IMS-Zubehör

Accessori IMS e Prodotti di Consumo

Accesorios de IMS

ACCESSOIRES IMS™

IMS™ SMALL/LARGE PARTS BOX




IMS™ Large Parts Box
280119
IMS-1273

IMS™ Small Parts Box
659401
IMS-1271

Designed to keep small items such as rubber dam clamps and matrix band retainers organized during the sterilization process.


IMS™ ENDODONTIC STAND

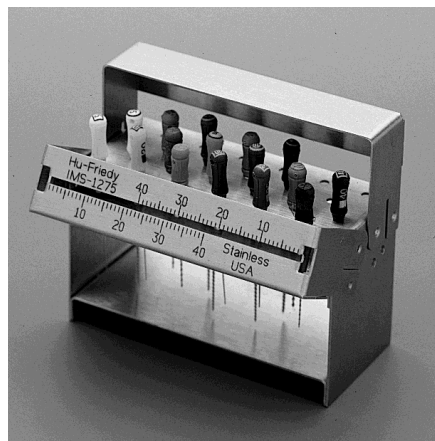

IMS™ Endodontic Stand
748366
IMS-1275

Hinged, double-lid stand holds up to 24 endodontic files and reamers. Lid is etched with a ruler to measure stop positions easily from either direction.


IMS™ BUR CUSHION

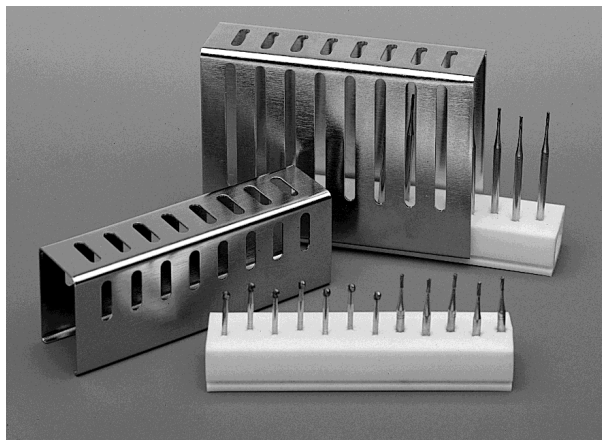

IMS™ Bur Cushion Short Lid
248231
IMS-1372S

IMS™ Bur Cushion Tall Lid
225126
IMS-1372T

A hollow, rectangular, silicone base with a stainless steel lid that is designed to hold 12 burs securely in place during cleaning and sterilization.

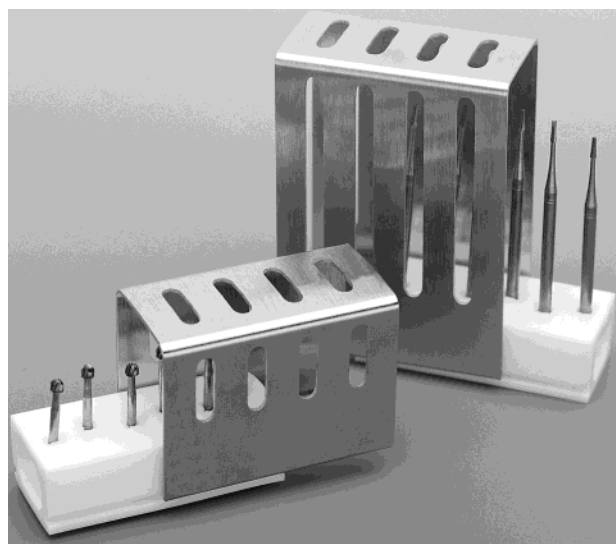

IMS™ Bur Cushion Tall Lid
588431
IMS-1372TH

IMS™ Bur Cushion Small Lid
997382
IMS-1372SH

Designed to hold 6 burs during cleaning and sterilization.

IMS Accessories

Accessoires IMS

Accessori IMS e Prodotti di Consumo

Accesorios de IMS

CLEANING AND CARE

IMS™ DAILY CLEAN



1lb/500g
makes 35
gallons/132l
173622
IMS-1220

5lb/2.27kg
makes 160
gallons/606l
198684
IMS-1218

Concentrated powder detergent specifically formulated for efficient ultrasonic cleaning or presoaking of cassette contents. Eliminates hand scrubbing.



A liquid cleaning concentrate on a Bi-enzyme base, ideal for the cleaning of instruments in an ultra-sonic bath or hand cleaning. Especially effective when cleaning instruments smeared with blood. The enzymes break down the proteins easily. Available in 1, 1 ltr pump bottles, in packs of 40 sachets of 10ml or in the large 3,8 ltr bulk container (not shown).

IMS™ PENETRATING OIL AND IMS INSTRUMENT LUBRICANT (NON-AEROSOL PUMP BOTTLE)




IMS™
Penetrating
Oil
931403
IPS

Use the penetrating oil to remove rust spots and corrosion.



IMS™
Instrument
Lubricant
136123
ILS

Makes for consistently smooth hinges of forceps, rongeurs, needle holders, hemostats etc. Apply after cleaning, prior to sterilization; it will not dry out during sterilization. Prevents corrosion and other deposits ensuring smooth operation.

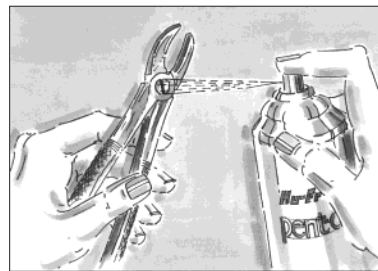

40 sachets
of 10ml
310565
IMS-1222


950ml
Pump-bottle
658944
IMS-1224

3,8l
bulk
container
336173
IMS-1226

- Cleaning and care
- Nettoyage et maintenance
- Pulizia e Cura
- Limpieza y Cuidado

STERILIZATION

IMS™ AUTOCLAVE WRAP

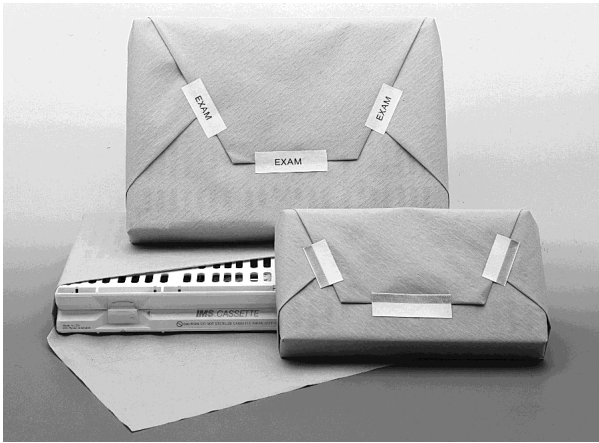

610 x 610mm
(24"x24") for
large cassettes
500 sheets
144401
IMS-1210

380 x 380mm
(15"x15") for
small cassettes
1000 sheets
325163
IMS-1211

500 x 500mm
(20"x20") for
1/2 DIN
Cassettes
152684
IMS-1212

Single-use wrap for sterile protection. When the cassette is unwrapped, the wrap provides a sterile environment for the cassette to be placed on. In addition, the wrap can be used to hygienically dispose of used cotton wool rolls and gauze pads. A wrapped cassette remains sterile up to 30 days.

IMS™ MULTIPLE TAPE DISPENSER



IMS™
Multiple Tape
Dispenser
164887
IMS-1270

Hold up to 9 rolls of IMS Monitor Tape. Rubber suction cups on the base secure dispenser to a table or counter. May also be wall mounted.

IMS™ STEAM STERILIZATION INDICATORS



IMS™ Steam
Sterilization
Indicators
173894
IMS-1240

Indicator strip verifies that steam sterilization parameters of time and temperature have been met. Use one in each cassette and retain for office records. 250 per box.

IMS™ MONITOR TAPE

Strong tape (3/4 x 60 yds, 2 cm x 55 m) for identifying cassettes (date, instrument group etc.) and visually differentiating processed from non-processed cassettes. Tape is available with specific procedural set-ups, color-coded or blank (see list below).

Blank	332755	IMS-1250
Amalgam	474657	IMS-1251
C&B	252395	IMS-1252
Cement	739634	IMS-1253
Composite	824275	IMS-1254
Endo	894176	IMS-1255
Exam	191221	IMS-1256
Extraction	345687	IMS-1257
Perio	164534	IMS-1258
Prophy	231215	IMS-1259
Pedo	343692	IMS-1260
Yellow	169455	IMS-1265
Red	979392	IMS-1266
Green	179418	IMS-1267
Blue	367249	IMS-1268
Black	568619	IMS-1269

- _____ Sterilization
- _____ Sterilisation
- _____ Sterilizzazione
- _____ Esterilización de IMS

ORGANIZATION

SMALL AND LARGE IMS COLOR-CODED INSTRUMENT RINGS



Silicone rubber rings for identifying and organizing instruments. Useful for maintaining proper instrument order in cassette and identifying instruments by procedure. Will withstand all heat sterilization methods.

Mini			Maxi	
Assorted Colors	845973	IMS-1280	545791	IMS-1380
Grey	147796	IMS-1281	--	
White	811971	IMS-1282	261084	IMS-1382
Orange	147338	IMS-1283	--	
Brown	886483	IMS-1284	--	
Yellow	949132	IMS-1285	164089	IMS-1385
Red	209589	IMS-1286	679259	IMS-1386
Green	306061	IMS-1287	533521	IMS-1387
Blue	153958	IMS-1288	234818	IMS-1388
Black	182753	IMS-1289	789615	IMS-1389

Small* color-coded rings, all one color, 100 per package.

Small* color-coded rings, assorted colors, 120 per package.

*Small = standard handle size.

Large* color-coded rings, all one color, 50 per package.

Large* color-coded rings, assorted colors, 60 per package.

*Large = e.g. used for forceps.

IMS™ CASSETTE I.D. LABELS AND I.D.



IMS™
Cassette
I.D. Labels
946532
IMS-1234

Specially designed polyester pressure-sensitive labels for permanent cassette identification, 50 per package.

IMS™
I.D.
Marker
245321
IMS-1235

Permanent indelible ink marker is ideal for identifying procedural set-ups. For use on I.D. Labels and Monitor Tape.

- Organization
- Organisation
- Organizzazione
- Organización de IMS

MAINTENANCE TIP

Be sure to date every cassette that is processed by writing on the monitor tape. This allows staff to easily rotate procedural set-up inventory.